



Electric Breast Pump

Hi bebe Plus

 **bistos**



Efficiency

Double Breast Pump

Don't let the size fool you. Hi-bebe plus is small and light enough to carry anywhere but doesn't compromise its performance. Dual expression and high pressure makes the expression in short time.

Comfort

Fine Suction and Cycle Control


10 step of suction and 3 step of cycle control per each mode helps expression at most comfort level. Expression doesn't mean to discomfort

More quiet than whispering sound

Quiet operating sound doesn't disturb mother and baby. Quiet sound brings peace to sleeping baby and mother



Hygiene

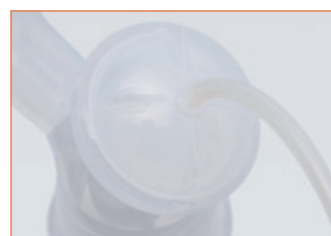

Closed pumping system

Generated air from the pump blocked by Diaphragm preventing infection. Also, this structure prevent backflow of milk into the tube and the pump.



BPA Free / Latex Free

All the parts contact the milk is made of BPA And Latex free material.

Portability

Rechargeable battery

Don't need to bring power adaptor for expression up to 150 min. Feel the freedom of electric cable.

Light Weight and Small Size

Babies are tiny but their necessities are not. Small size and light weight helps mothers carry the pumps anywhere with baby care product.



Convenience

Smart memory mode

Once you set your frequent use of pressure level and the time of each mode, Turning on power is only thing to do. Mothers can concentrate on other work. Up to 40 sequence can be saved.

Easy to use

Assembling & dis-assembling, clean, dry is very easy and only takes short time.



BISTOS Technical Specification

Model	BT - 100	
Order code	100S	100D
Funnel & Bottle set	1 set	2 set
Funnel Size	24mm (Option: 21 / 27mm)	
Materials	BPA Free / Latex Free	
Display	3" Mono LCD	
Mode	Expression, Massage, Memory	
Vacuum Level	-20 to -270mmHg	
Vacuum Control	10 Step	
Expression Cycle	24 to 70 Cycle/Min (Expression)	
Cycle Control	3 Step	
Auto shut off	30 Min.	
Pumping system	Close system (Backflow prevention)	
Battery Type (Capacity)	Li-ion (900mA)	
Battery Run Time	150 Min.	
Battery Charging Time	90 Min.	
Optional Accessories	Carrying Bag, Car Charging Cable, Funnel & Bottle set	
Warranty	1 year (Main unit)	
Dimension (HxWxD)	132 x 95 x 49mm	
Weight	260g	
Certificates	FDA, CE, ISO13485, GMP	
Origin	Rep. of Korea	



Carrying Bag



Single pump



Funnel 21 / 24 / 27mm

BIO SIGNAL TOTAL SOLUTION

Bistos Co., Ltd. (Headquarter)

7th Fl., A Bldg., Woolim Lions Valley 5-cha, 302,
Galmachi-ro, Jungwon-gu, Seongnam-si,
Gyeonggi-do, Korea
Tel. +82 (31) 750 0340
Fax. +82 (31) 750 0344

bistos@bistos.co.kr
www.bistos.co.kr

Bistos America Inc. (US Branch)

11417 W Bernardo Ct Ste D, San Diego, CA 92117, USA
Tel. +1-868-433-6689
Fax. +1-888-391-5153