

FUNCTIONS



Backrest adjustment : 0 - 68°



Thigh rest adjustment : 0 - 27°



Foot rest mechanical adjustment : 0 - 22°



Vascular position




Fowler position




Trendelenburg position : 0 - 12°

OPTIONS

GROUP 1 / CASTOR OPTIONS

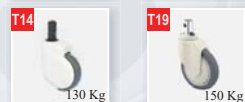


GROUP 2 / SIDE RAIL OPTIONS

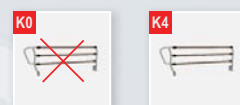

GROUP 3 / ELECTRICAL ACCESSORY OPTIONS



GROUP 4 / MATTRESS PLATFORM OPTIONS

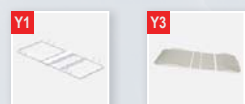

GROUP 5 / OTHER ACCESSORY OPTIONS



DESCRIPTION OF OPTION CODES

T1 TWIN CASTOR TWO WITH BRAKE (120 MM)	E3 UNDER BED LIGHTING SYSTEM
T16 CASTOR TWO WITH BRAKE (125 MM)	Y1 REMOVABLE 4 PARTED MESH PLATFORM.
T4 METAL CASTOR TWO WITH BRAKE (125 MM)	Y3 REMOVABLE DOUBLE WALLED PE BED PLATFORM
T14 CASTOR CENTRAL & DIRECTIONAL LOCK (125 MM)	W0 WITHOUT OTHER ACCESSORIES
T19 CASTOR CENTRAL & DIRECTIONAL LOCK (150 MM)	W1 IV STAND : ONE HEIGHT ADJUSTMENT SYSTEM
K0 WITHOUT SIDE RAIL	W2 LIFTING POLE : BEARING CAPACITY 75 Kg
K4 WITH ALUMINIUM SIDE RAIL	W4 TELEVISION STAND FOR 17" LCD TV
E0 WITHOUT ELECTRICAL ACCESSORIES	W5 TRACTION APPARATUS
E1 BATTERY RECHARGEABLE LEAD GEL TYPE. OUTPUT VOLTAGE: 24 V DC. CAPACITY: 1,2 Ah.	

T1, T16, T4, T14, T19, K0, K4, E0, E1 Must select only one of them.

E3, Y1, Y3, W1, W2, W4, W5 Do not select the others if you select without accessory. If you select the others you can select all or 1-2-3 ... of them.

GENERAL SPECIFICATIONS

- Easy to clean four-piece bed platform.
- Head and foot movements with the help of an electrical motor.
- Trendelenburg movement by the help of shock absorber mechanism at foot side.
- Ergonomic, durable and easily cleanable aluminum rails system with maximum safety designed in accordance with the standards of EN 60601-2-52.
- Head and foot platforms with automatic sliding system.
- Easily removable and installable, hygienic, high resistant head board system designed of compact laminate and aluminum.
- Shock resistance plastic bumper wheels.
- Metal components dyed with anti-corrosion electrostatic powder paint after being treated with phosphate.
- Optional laminate color choice.^(RLX)

TECHNICAL SPECIFICATIONS

■ External dimension	: 102 x 222,5 x 49,5 cm
■ Width outside side rails	: 98 cm
■ Mattres platform	: 89 x 195 cm
■ Recommended mattress dimension	: 85 x 194 x 12 cm
■ Autoregression (Backrest / Legrest)	: 8 cm / 4 cm
■ Safe Operating Load (Inc. Mattres and Accessories)	: 230 kg
■ Weight	: ~ 101 kg
■ Volume	: 1,23 m ³

ELECTRICAL SPECIFICATIONS

■ Power Input	: 230 V AC 50 Hz Max 1 A
■ Power Output	: 24 V DC
■ Electrical Current Protection	: Class 1, Type B
■ Liquid and Dust Protection	: IP54 – IP66



* RL : Laminate color options page 84.